**intensive massage with red light**  
**massage intensif avec lumière rouge**  
**masaje intenso con luz roja**  
**massagem intensiva**  
**com calor por luz vermelha**

prorelax®

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# intensive massage with red light

## Table of Contents

Introduction	3
Features	3
Safety Instructions	4-5
Precautions	6
Product Description	7
Instructions for use	8-9
Technical Data / Maintenance	10
Warranty Provisions	11
Contact	42

**prorelax**<sup>®</sup>

Dear **prorelax**<sup>®</sup> Customer,

thank you for purchasing our **prorelax**<sup>®</sup> INTENSE BODY MASSAGER.

We want you to be completely satisfied with this product's quality and performance. Many people, especially those who play golf or tennis, are experiencing pain and muscular tension, which are definitely taking the fun out of the sports.

Our **prorelax**<sup>®</sup> INTENSE BODY MASSAGER should help you to overcome those difficulties. In order to take full advantage of this product's features and benefits, we ask you to carefully read these instructions.

## FEATURES

The **prorelax**<sup>®</sup> INTENSE BODY MASSAGER is a high quality product with innovative features different from older percussion massagers:

- Light weight with extra long ring type detachable handle for easy reach of all parts of the body.
- Large heat surface for soothing of muscles.
- Variable massage strength selection.
- Incorporated hand strap for palm gripping of massager.
- Two attachments for special massage.

When using an electrical appliance, basic precautions should always be followed.

## **Read all instructions before using this appliance!**

### **To reduce the risk of electric shock:**

1. Always disconnect the appliance from the electrical outlet immediately after use and before cleaning.
2. DO NOT grab for an appliance which has dropped into water. Unplug it immediately.
3. DO NOT use the massager while in a bath tub or in a shower.
4. DO NOT place or store the appliance where it can fall or be pulled into a tub or sink.
5. **DO NOT place the massager in or drop into water or other liquid.**

### **To prevent the risk of burns, fire, electric shock or injury, please make sure that:**

1. the appliance will never be left unattended when plugged in. Unplug from outlet when not in use and before adding or removing attachments.
2. Do not operate it under a blanket or pillow. Excessive heating can occur and cause fire, electrical shock or injury.
3. Close supervision is necessary when this appliance is used on, or near children, or disabled persons.
4. Use this appliance only for its intended purpose as described in this manual.  
Do not use attachments not provided with the unit.

5. Never operate this appliance if it has a damaged cord or plug, or if not working properly, or if it has been dropped and damaged, or dropped into water. Return it to the Service Center for examination and repair.
6. Do not hold this appliance by the electrical cord or use the cord as a handle.
7. Keep cord away from heated surfaces.
8. Never drop or insert anything into any opening.
9. To disconnect, turn all controls to the "OFF" position, then remove the plug from the wall outlet.
10. Do not use the massager outdoors.
11. Do not operate where aerosol (spray) products are used or where oxygen is being administered.
12. Do not step on the massager.
- 13. Never use the massager with the air openings blocked. Keep the air openings free of lint, hair, and the like.**
14. Use heated surface carefully. It may cause serious burns. Do not use on **areas** of sensitive skin or if you have poor circulation. Never use heat on children or disabled persons.
15. **If your massager does not work properly anymore, please do not try to fix it by yourself. Only an authorized dealer should do this in order to make sure the device is being serviced correctly and will be safe for use again. Unauthorized attempts to repair the product will void your warranty.**

## Precautions

The **prorelax**® INTENSE BODY MASSAGER is a state-of-the-art massaging device that needs to be operated in compliance with this manual. There are some restrictions where and when not to use it:

The red light massager should not be used in the following cases: Tuberculosis, benign and malignant tumors, hemorrhages, inflammations of the skin, phlebitis and thrombosis, open and fresh wounds, bruises, broken skin or varicose veins, undetermined calf pain, or joint pain.


In case of doubts, please ask you physician before using it.

Do not use if you are experiencing unexplained pain or swollen muscles, or following a serious muscular injury before consulting your physician.

Discontinue use of massager if pain increases or swelling occurs after application. If you are pregnant, please consult your physician in order to find out whether a massage is good for you or not.

## PRODUCT DESCRIPTION

- 1) Massage/massage + red light heat/off control switch
- 2) Variable massage strength control knob
- 3) Velcro hand grip strap
- 4) Red light heat massage head
- 5) Ring-type handle (detachable)
- 6) Handle release buttons
- 7) Massage body cover
- 8) Attachment for the use of lotion or cream
- 9) Accunodes attachment to intensify the massage and the blood circulation

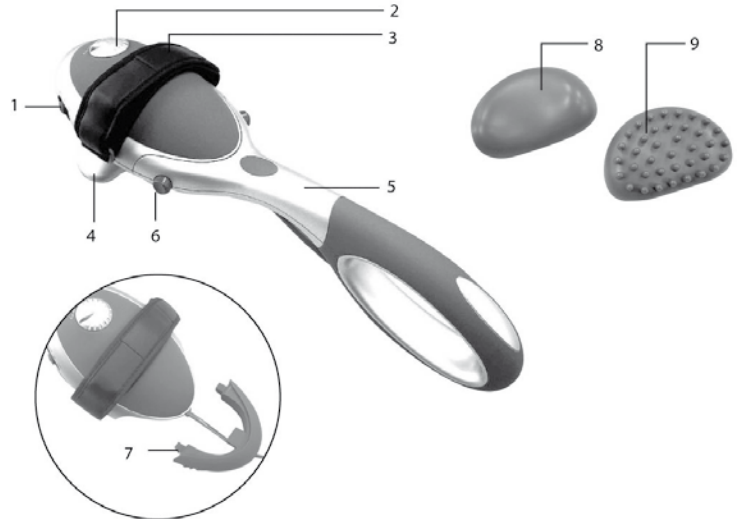---

## Instruction for Use

- Make sure that the massager is set to **“0, off”** position before inserting the plug into a 120V household electrical outlet.
- To turn unit on, simply slide the power switch to **“M & M+H”** position.
  1. M = Massage only
  2. M+H = Massage red light heat
- You can attach the two massage attachments for different massage applications. For a strong and hard massage, use the fixed round massage head.

The acupressure massage attachment will provide a larger massage area, and can be used with crèmes and lotion, while the acunodes attachment will provide a distinct multi-point massage that is more intense for increased blood flow.

We recommend not massaging one area for more than 5 minutes. The actual massage time depends on your own feeling.

- We recommend not using the massager for more then 20 minutes continuously as it has to cool down between applications.
- When massage is completed turn off the appliance by sliding to the **“0”** position and unplug and store it in a dry place.
- This appliance is for household use only.
- For selecting the massage strength, rotate the **“variable strength control knob”** with your finger slowly from **“-”** to **“+”** until you get to the desired strength.

### Detachable ring-type handle & palm-type feature

This red light palm percussion massager does have a detachable handle and palm-type handling feature.

To use the **“palm type handling”**, simply press the two “handle release buttons” and pull out the handle. Use the “body cover” to close the opening (see photo) and adjust the Velcro hand strap. Now you can use the massager in a convenient way for the massage of various parts of the body, especially the arms, thighs and calves. For massaging hard to reach body areas reattach the ring-type handle back by reversing the above procedure.

## TECHNICAL DATA

Name and model	HAND MASSAGER WITH RED LIGHT PSA-2365
Power supply (Europe)	230-240V~50Hz
Power input (Europe)	22 W
Power supply (US)	120 V ~ 60Hz
Power input (US)	14 W
Protection class	II
Short-term operation	20 minutes
Ambient conditions	Only for dry application and private use
Dimensions (L x W x H)	40 x 11.5 x 12 cm

## MAINTENANCE

### TO CLEAN:

Be sure to unplug the unit and allow it to cool before cleaning. Use a soft, slightly damp cloth to wipe. Never allow water or any liquids to come into contact with the unit.

### TO STORE:

Unplug the appliance from the outlet and allow it to cool before storing it in its box or a clean, dry place. Do not wrap the power cord around the unit. This can cause the cord to wear and break. Do not hang the unit by the power cord.

## WARRANTY PROVISIONS

We guarantee completeness and freedom from defects at the time the sale is made.

Repairs may be undertaken exclusively by importers, service centers or dealers who have been authorized to do so by PLI.

The warranty period shall be two years from the date of purchase subject to the requirement that use is made as described in the usage directives. Parts which are subject to normal wear-and-tear (e.g., batteries) are not covered by the warranty.

Defects discovered within the warranty period must be reported directly to the authorized dealer or the service center so that the cause can be determined. If it is not possible to repair the defects, then the product shall be exchanged.

The timeframe required for making repairs shall be dependent upon the type of defect.

The Buyer shall assume the costs for the delivery and pick-up of the device.

The warranty to be provided by PLI shall be limited to the amount of the purchase price. The Buyer shall be entitled to no more extensive claims.

**In the following cases, the warranty claim shall not apply:**

- a) No receipt can be presented.
- b) The Buyer or a non-authorized person has attempted to repair the device.
- c) The defect was caused through improper usage upon the part of the Buyer, particularly through non-adherence to the usage directives.
- d) In the event of obvious damage caused by fluids, by the device being dropped or
- e) through blunt force exerted during the usage of the product.

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PRORELAX INTERNATIONAL LLC
PROMEDICS HEALTHCARE GmbH